



DC150 DC250 DC350 Blanket Warming Cabinet

37-82°C (98-180°F)
temperature range

WarmRight® technology
uses multiple zone heating
that monitors surfaces for
safe temperature ranges.

Efficient, balanced control
allows for low energy
consumption and
minimal heat loss.

Interior lighting casts
a comforting blue glow
with two (2) different
intensity settings.

Epoxy-coated
interior surfaces

Available in three (3)
different cabinet sizes



DC350



DC250



DC150

 **ENTHERMICS**
Medical Systems

www.enthermics.com

*Experience the Enthermics difference:
Quality beyond comfort.*

DC150, DC250, DC350 Blanket Warming Cabinet

Specifications:


- Single-chamber warming cabinet
- White epoxy-coated steel exterior casing and interior insert
- Single pane, energy efficient, e-coated glass window in door allows for inventory observation
- Easy, hands-free, push-button door release
- Door is fully gasketed and hinged on the right side of the unit
- WarmRight® technology utilizes multiple zone heating that monitors surfaces for safe temperature ranges
- DC350 includes a heated center shelf
- Furnished with four (4) 3/4" (19mm) non-skid rubber feet
- Optional casters available

Control:

- Electronic control
- Efficient, balanced heating system allows for low energy consumption and minimal heat loss
- Adjustable temperature range of 37 - 82°C (98 - 180°F)
- Operates in Celsius or Fahrenheit
- Four digit L.E.D. display
- On/Off key
- Up and down arrow keys
- Actual temperature key
- Interior light key
- Built-in speaker for audible feedback
- Built-in button lock-out feature on control

Additional features:

- L.E.D. interior lighting casts a comforting blue glow with two (2) different intensity settings, and off mode.
- In the event of a power failure the cabinet will remember its programming and begin to operate as before when power is restored.
- A warming shut-off system, separate from the electronic control, prevents overheating.
- Optional timer allows the user to program the control to automatically turn on and off at a selected time.

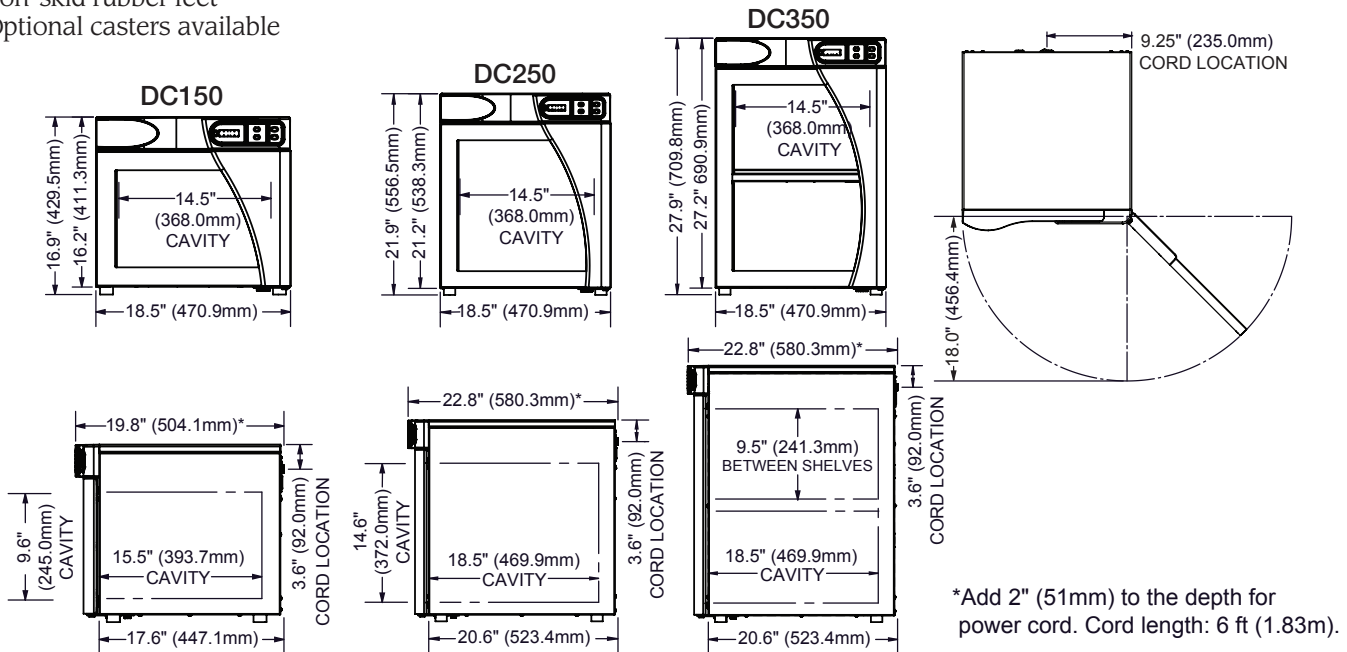

Model	DC150	DC250	DC350
Interior Capacity	1.5 ft ³	2.5 ft ³	3.5 ft ³
Net Weight	56 lbs (25 kg)	67 lbs (30 kg)	91 lbs (41 kg)
Ship Weight	58 lbs (26 kg)	135 lbs (61 kg)	167 lbs (76 kg)
Electrical			
120 V.A.C. 50/60 Hz, 1ph	5.0 Amps, 0.6 kW NEMA 5-15P 15A - 125V Plug Hospital Grade	5.0 Amps, 0.6 kW NEMA 5-15P 15A - 125V Plug Hospital Grade	6.7 Amps, 0.8 kW NEMA 5-15P 15A - 125V Plug Hospital Grade
230 V.A.C. 50/60 Hz, 1ph	2.6 Amps, 0.6 kW BS 1363 Plug (UK only) CEE 7/7, 220-230V Plug	2.6 Amps, 0.6 kW BS 1363 Plug (UK only) CEE 7/7, 220-230V Plug	3.4 Amps, 0.8 kW BS 1363 Plug (UK only) CEE 7/7, 220-230V Plug

Options & Accessories	
Casters, 3" (25mm)	5012693
Timer	Available with original order only



www.enthermics.com
An ISO 13485:2003 certified company

W164 N9221 Water Street • PO Box 443
Menomonee Falls, WI 53052-0443
Phone 262.251.8356 • 800.862.9276
Fax 262.251.7067
Email generalinfo@enthermics.com