



EC1540

Blanket Warming Cabinet

15.4 ft³ (0.44m³) interior capacity

Glass doors



Shown with optional
full-perimeter rubber bumper

37-94°C (98-200°F) range

Deluxe electronic control

24-hour timer
(single programmable "off" period)

Plastic-coated interior surfaces

 **ENTHERMICS**
Medical Systems

www.entermics.com

*Experience the Entermics difference:
Quality beyond comfort.*

EC1540 Blanket Warming Cabinet

Specifications:


Dual-cabinet design with 20 gauge stainless steel exterior casing and doors. Window doors allow observation of inventory with the doors closed. The cabinet is warmed using a low-heat-density electrothermal cable array mounted against the walls, eliminating the need for a heat circulating fan. Each compartment is equipped with one (1) white, epoxy-coated blanket support assembly and one (1) center shelf. The cabinet is furnished with one (1) set of 5" (127mm) heavy duty casters, two with locking brake.

Special Information:

- Right-hand door swing is standard. Left-hand door is available as a special order.
- Pass-through design is available.

Control:

Each compartment is controlled by one (1) electronic control which consists of a four digit L.E.D. display, on/off key, heat indicator L.E.D., increase and decrease keys, integrated lock feature and a series of prompt sequence indicators. Each control has an adjustable temperature range of 37 - 94°C (98 - 200°F). An alarm will sound if temperatures exceed 10°F (6°C) above the set-point temperature. Controls can easily be set to operate in Celsius or Fahrenheit. Timer key allows the user to program the control to automatically turn on and turn off once at a selected time during a 24 hour period. In the event of a power failure, the cabinet will remember its programming and begin to operate as before when power is restored. A warming shut-off system, separate from the electronic control, prevents overheating.

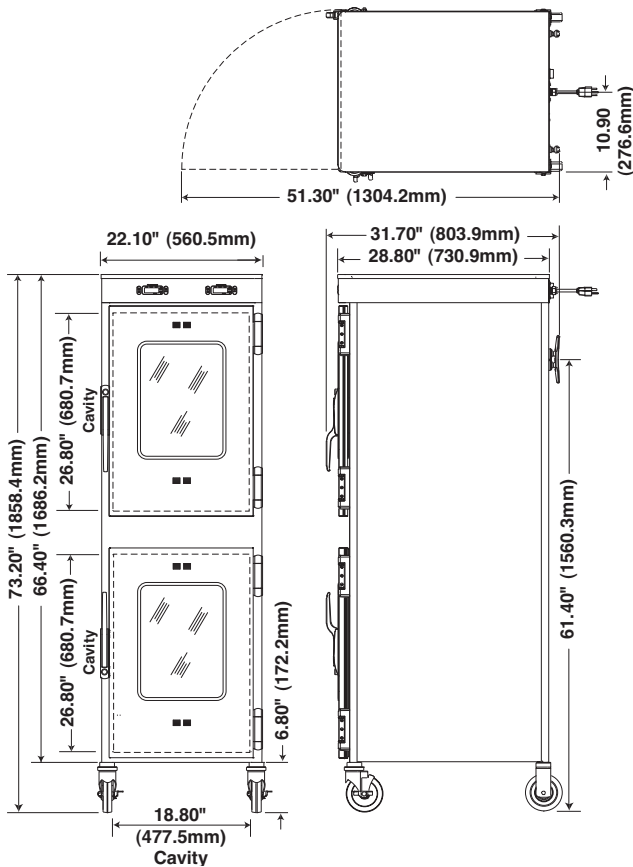

Capacity and Weights	
Interior Capacity per Compartment	7.7 cubic feet (c. 0.1m ³)
Net Weight	412 lb (187 kg)
Ship Weight	515 lb (234 kg) est.
Electrical	
120 V.A.C. 50/60 Hz, 1ph 1920 Watts, 16.0 Amps	NEMA 5-20P 20A - 125V Plug Hospital Grade
230 V.A.C. 50/60 Hz, 1ph 1700 Watts, 7.4 Amps	BS 1363 Plug (UK only)
	CEE 7/7 220-230V Plug
Options and Accessories	
Cylinder Lock for Door Handle	LK-22567
Full-Perimeter Rubber Bumper Assembly	5001763
Leg Kit, 6" (152mm)	44093
Window Door Assembly, LH, ea.	5001255
Solid Door Assembly, RH, ea.	5005808
Solid Door Assembly, LH, ea.	5006023
Built-in Trim Kit	5003168



**ENTHERMICS
Medical Systems**
www.enthermics.com

W164 N9221 Water Street • PO Box 443
Menomonee Falls, WI 53052-0443
Phone 262.251.8356 • 800.862.9276
Fax 262.251.7067
Email generalinfo@enthermics.com