



EC2060

Blanket Warming Cabinet

20.6 ft³ (0.58m³) interior capacity

Glass doors



37-94°C (98-200°F) range

Deluxe electronic control

24-hour timer
(single programmable "off" period)

Plastic-coated interior surfaces

 **ENTHERMICS**
Medical Systems

www.enthermics.com

*Experience the Entermics difference:
Quality beyond comfort.*

EC2060 Blanket Warming Cabinet

Specifications:


The single-chamber, dual door warming cabinet is constructed with 20 gauge stainless steel exterior casing and double door design. Window doors allow observation of inventory with the doors closed. The cabinet is warmed using a low-heat-density electrothermal cable array mounted against the walls, eliminating the need for a heat circulating fan. The warming cabinet is equipped with one (1) white, epoxy-coated blanket support assembly and one (1) center shelf per door opening. The cabinet is furnished with one (1) set of 5" (127mm) heavy duty casters, two with locking brake.




Special Information:

- Right-hand door swing is standard.
Left-hand door is available as a special order.
- Pass-through design is available.

Control:

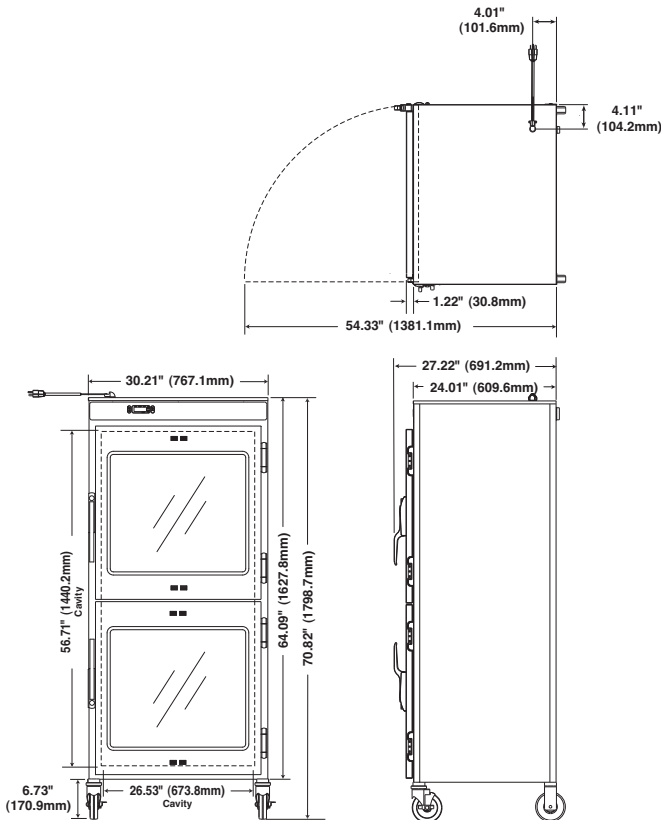The chamber is controlled by one (1) electronic control which consists of a four digit L.E.D. display, on/off key, heat indicator L.E.D., increase and decrease keys, built-in lock-out feature and a series of prompt sequence indicators. The control can easily be set to operate in Fahrenheit or Celsius and has an adjustable temperature range of 37 - 94°C (98 - 200°F). Timer key allows the user to program the control to automatically turn on and turn off once at a selected time during a 24 hour period. In the event of a power failure, the cabinet will remember its programming and begin to operate as before when power is restored. A warming shut-off system, separate from the electronic control, is an additional safety feature

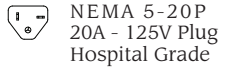

Capacity and Weights	
Interior Capacity per Compartment	20.6 cubic feet (c. 0.58m ³)
Net Weight	356 lb (161 kg)
Ship Weight	485 lb (220 kg) est.
Electrical	
120 V.A.C. 50/60 Hz, 1ph 1700 Watts, 14.2 Amps	  
Options and Accessories	
Leg Kit, 6" (152mm)	44093
Built-in Trim Kit	5009266



**ENTHERMICS
Medical Systems**
www.enthermics.com

W164 N9221 Water Street • PO Box 443
Menomonee Falls, WI 53052-0443
Phone 262.251.8356 • 800.862.9276
Fax 262.251.7067
Email generalinfo@enthermics.com